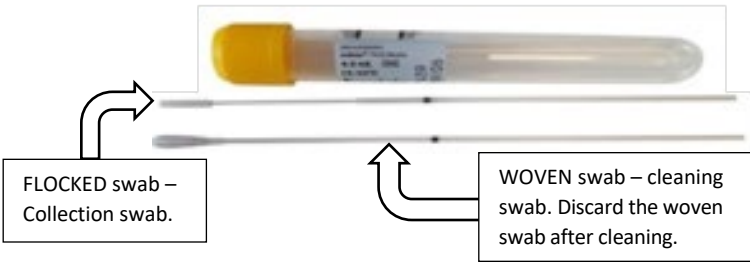


UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
HER-2 FISH (Primary/Metastatic Breast only)	FFPE Only	
- Factor V Leiden Mutation - Prothrombin (Factor II) DNA Mutation	K2 EDTA	
- Hereditary Hemochromatosis (HFE) - Cystic Fibrosis DNA Mutations	EDTA	
Viral Loads (Quantitative PCR) - HIV-1 - CMV - HBV - HCV - EBV - BK -Blood	Plasma from EDTA	
BK Viral Load – Urine	Cobas PCR Media – Urine Collection Kit	
- Chlamydia Trachomatis - Neisseria Gonorrhoeae (CT/NG) – URINE	Cobas PCR Media – Urine Collection Kit	

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

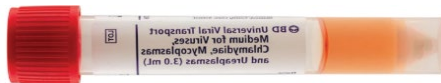
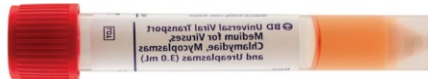
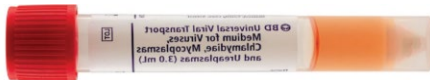
TEST	SPECIMEN COLLECTION CONTAINER
Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) FEMALE GENITAL – ENDOCERVICAL COLLECTION	Cobas PCR Media Dual Swab Kit ONLY (refer to kit instructions) 
Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) FEMALE GENITAL – VAGINAL COLLECTION (Clinician and Self-Collections)	Cobas PCR Media Uni or Dual Swab Kit (Woven swab ONLY) 
Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) ANORECTAL COLLECTION (Clinician and Self-Collections)	Cobas PCR Media Uni or Dual Swab Kit (Woven swab ONLY) 
Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) THROAT COLLECTION (Clinician and Self-Collections)	Cobas PCR Media Uni or Dual Swab Kit (Woven swab ONLY) 

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
Human Papillomavirus (HPV) DNA CERVICAL COLLECTION (Clinician Collections)	Hologic ThinPrep PAP Test Vial <i>Cervical collections using PAP collection devices (spatula brush, or broom) into Hologic ThinPrep PAP Test vial. Refer to Lab Test Directory for detailed instructions.</i> <i>Note: ThinPrep PreservCyt Solution (NON-GYN) vial is unacceptable for HPV DNA testing.</i>	 
Human Papillomavirus (HPV) DNA VAGINAL COLLECTION (Patient Self-Collections only) <i>Note: Clinician-collected vaginal specimens using Copan FLOQSwabs are unacceptable and will be rejected.</i> <i>Note: Clinician-collected vaginal specimens using PAP collection devices (spatula, brush, or broom) will be forwarded to send out reference lab when necessary.</i>	Hologic ThinPrep PAP Test Vial <i>Vaginal self-collections using Copan FLOQSwab 552C into Hologic ThinPrep PAP Test vial. Refer to Lab Test Directory for detailed instructions.</i> <i>Note: ThinPrep PreservCyt Solution (NON-GYN) vial is unacceptable for HPV DNA testing.</i>	Copan FLOQSwab 552C 1. Lawson Stores 103 Inventory Template #198951 (box of 100) 2. Lawson Catalog #198678 (individual swabs or case of 1000)  

TO ENSURE ACCURACY OF MOST UP TO DATE INFORMATION, PLEASE DO NOT PRINT & DISTRIBUTE

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
Respiratory Viral Panel	BD UTM	
	M4	
	M4 RT	 
Covid-2019 RNA Qual	BD UTM	
	M4 RT	 
Covid-2019 + Flu A/B	BD UTM	

TO ENSURE ACCURACY OF MOST UP TO DATE INFORMATION, PLEASE DO NOT PRINT & DISTRIBUTE

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

For additional details regarding test requisitioning, specimen collection and transport requirements, results, and interpretation, search individual tests in our Lab Test Directory at

<https://www.testmenu.com/ucdavis>

Any further questions should be routed to Client Services Department:

Phone: (916) 734-7373 Fax: (916) 737-1001