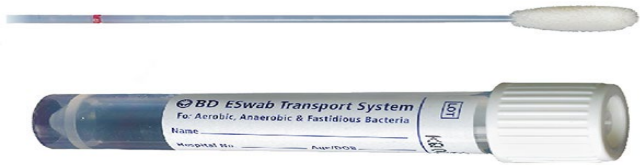
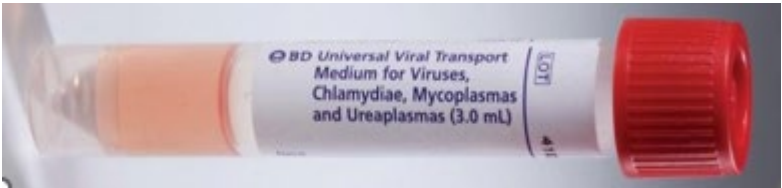


MICROBIOLOGY COLLECTION GUIDE

Bacterial Culture Swab (Dual Red) PS# 15347 	Group A Strep Rapid Antigen and reflex culture- <u>NO ALTERNATIVE</u> Also acceptable, but NOT preferred for : • Aerobic cultures-wound, genital, nose, etc. • MRSA surveillance
E-Swab PS# 72179 	Preferred swab for routine cultures: • Anaerobic culture • Aerobic cultures-wound, genital, nose, etc. • MRSA surveillance • Wet Preps • Yeast culture-throat, mouth, genital
Universal Viral Transort PS# 69846 	Required swab for testing, <u>NO ALTERNATIVE</u> • Bordetella pertussis PCR • HSV/VZV Combo (Herpes/Varicella) • Covid, Influenza, or RSV • RVP-Respiratory Viral Panel • Mini-RVP (Covid, Flu A+B, RSV) • NOTE: UVT cannot be used for bacterial routine cultures

MICROBIOLOGY COLLECTION GUIDE

Aptima Unisex Swab (foil cap)

PS# 72215



Available Tests

- Chlamydia trachomatis
- Neisseria gonorrhea
- Trichomonas vaginalis
- Mycoplasma genitalium
- Vaginitis Panel

Aptima Multitest Swab (foil cap)

PS# 72217



Available Tests

- Chlamydia trachomatis
- Neisseria gonorrhea
- Trichomonas vaginalis
- Mycoplasma genitalium
- Vaginitis Panel

Urine

PS# 15694



Available Tests

- Chlamydia trachomatis
- Neisseria gonorrhea
- Trichomonas vaginalis
- Mycoplasma genitalium
- Urine culture