









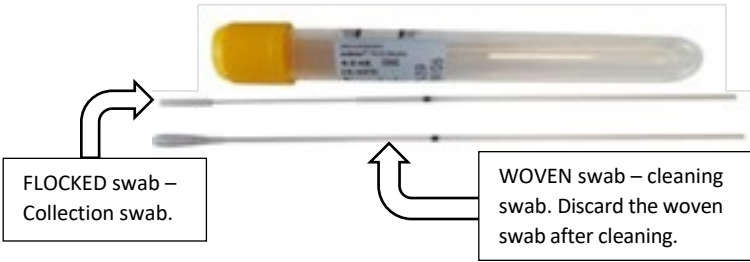



TO ENSURE ACCURACY OF MOST UP TO DATE INFORMATION, PLEASE DO NOT PRINT & DISTRIBUTE

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
<ul style="list-style-type: none"> • CLL FISH Panel 	Sodium Heparin	
	EDTA	
<ul style="list-style-type: none"> • HER-2 FISH (Breast & Gastric) 	FFPE Only	
<ul style="list-style-type: none"> • Factor V Leiden Mutation • Prothrombin (Factor II) DNA Mutation 	K2 EDTA	
<ul style="list-style-type: none"> • Hereditary Hemochromatosis (HFE) • Cystic Fibrosis DNA Mutations 	EDTA	
<p>Viral Loads (Quantitative PCR)</p> <ul style="list-style-type: none"> • HIV-1 • CMV • HBV • HCV • EBV • BK -Blood 	Plasma from EDTA	
<ul style="list-style-type: none"> • BK Viral Load – Urine 	Cobas PCR Media – Urine Collection Kit	
<ul style="list-style-type: none"> • Chlamydia Trachomatis • Neisseria Gonorrhoeae (CT/NG) – <p>URINE</p>	Cobas PCR Media – Urine Collection Kit	


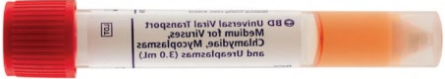


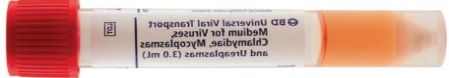
TO ENSURE ACCURACY OF MOST UP TO DATE INFORMATION, PLEASE DO NOT PRINT & DISTRIBUTE

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
<ul style="list-style-type: none"> Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) <p>FEMALE GENITAL – ENDOCERVICAL COLLECTION</p>	<p>Cobas PCR Media Dual Swab Kit ONLY <i>(refer to kit instructions)</i></p>	 <p>FLOCKED swab – Collection swab.</p> <p>WOVEN swab – cleaning swab. Discard the woven swab after cleaning.</p>
<ul style="list-style-type: none"> Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) <p>FEMALE GENITAL – VAGINAL COLLECTION</p>	<p>Cobas PCR Media Uni or Dual Swab Kit (Woven swab ONLY)</p>	 <p>WOVEN swab – collection swab</p>
<ul style="list-style-type: none"> Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) <p>ANORECTAL COLLECTION</p>	<p>Cobas PCR Media Uni or Dual Swab Kit (Woven swab ONLY)</p>	 <p>WOVEN swab – collection swab</p>
<ul style="list-style-type: none"> Chlamydia Trachomatis / Neisseria Gonorrhoeae (CT/NG) <p>THROAT COLLECTION</p>	<p>Cobas PCR Media Uni or Dual Swab Kit (Woven swab ONLY)</p>	 <p>WOVEN swab – collection swab</p>


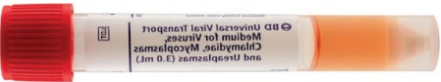


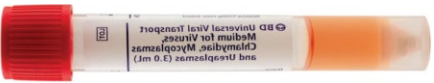
TO ENSURE ACCURACY OF MOST UP TO DATE INFORMATION, PLEASE DO NOT PRINT & DISTRIBUTE

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
<ul style="list-style-type: none"> Human Papillomavirus (HPV DNA) 	<p>Hologic ThinPrep PAP TEST</p>	<p><i>Important Notes:</i></p> <ol style="list-style-type: none"> Endocervical Brush and Spatula collections ONLY. Broom collection device is unacceptable for HPV DNA test on the Roche 4800 test system. ThinPrep PreservCyt Solution (NON-GYN) vial is unacceptable for HPV DNA test 
<ul style="list-style-type: none"> Respiratory Viral Panel 	<p>BD UTM</p>	
	<p>M4</p>	
	<p>M4 RT</p>	
<ul style="list-style-type: none"> HSV 1&2/VZV DNA, Swab 	<p>BD UTM</p>	

TO ENSURE ACCURACY OF MOST UP TO DATE INFORMATION, PLEASE DO NOT PRINT & DISTRIBUTE

UCDH MOLECULAR ACCEPTABLE SPECIMEN COLLECTION CONTAINERS

TEST	SPECIMEN COLLECTION CONTAINER	
<ul style="list-style-type: none"> • HSV 1&2/VZV DNA, Swab (<i>cont'd</i>) 	M4 RT	
<ul style="list-style-type: none"> • Bordetella Pertussis/Parapertussis DNA 	BD UTM	
	M4	
	M4 RT	
<ul style="list-style-type: none"> • HSV 1&2 DNA, CSF • Meningitis/Encephalitis Panel 	Sterile Tube (CSF from Lumbar Puncture)	
<ul style="list-style-type: none"> • Covid-19 RNA Qual • Covid-19 + Flu A/B 	BD UTM	
<p align="center">For additional details regarding test requisitioning, specimen collection and transport requirements, results, and interpretation, search individual tests in our Lab Test Directory at https://www.testmenu.com/ucdavis Any further questions should be routed to Client Services Department: Phone: (916) 734-7373 Fax: (916) 737-1001</p>		