



**STERILE CUP**

Acid-Fast Bacillus (AFB) Culture and AFB Stain (CXAFB)  
 ANAEROBIC CULTURE (CXANA)  
 ASPERGILLUS GALACTOMANNAN AG, BRONCH WASH (ASPGW)  
 BACTERIAL CULTURE AND GRAM STAIN (CXBS) –**COPAN eSWAB PREFERRED**  
 BACTERIAL CULTURE WITH ANAEROBE AND GRAM STAIN (CXBSN) –**COPAN eSWAB PREFERRED**  
 Bronchial Washing Culture, Quantitative (CXQBR)  
 Clostridioides (Clostridium) difficile Toxin PCR, Stool(CDPCI)  
 CHLAMYDIA/GC PCR-URINE (CTGCU)- **COBAS PCR URINE SAMPLE KIT PREFERRED**  
 FECAL LACTOFERRIN, STOOL WBC (LEUKV)  
 FUNGAL CULTURE (CXFUN)  
 FUNGAL CULTURE HAIR, NAIL, SKIN (CXFNS)  
 Giardia/Cryptosporidium Antigen, Stool (GCEIA)  
 Helicobacter pylori Antigen, Stool (HELO)  
 Histoplasma Antigen, Urine (HGU)  
 LEGIONELLA AG URINE (LEGAT)  
 Mycobacterium tuberculosis/Rifampin Resistance PCR (TBRIF)  
 PNEUMOCYSTIS JIROVECI PNEUMONIA (PJP)  
 Pneumonia Pathogen Nucleic Acid Detection Test (BALPC)  
 Tissue Culture, Quantitative (CXQNT)  
 QUANTITATIVE SMALL BOWEL CULTURE (CXQSB)  
 RESPIRATORY CULTURE W/GRAM STAIN (CXRES)  
 Sterile Body Fluid Culture with Gram Stain (CXFLD)  
 Sterile Body Fluid Culture with Anaerobe and Gram Stain (CXFLN)  
 TIP CULTURE (CXTIP)  
 TISSUE CULTURE W/GRAM STAIN, AEROBIC (CXTIS)  
 Tissue Culture with Anaerobe and Gram Stain (CXTIN)  
 URINE CULTURE (CXURN) – **BD VACUTAINER PREFERRED**  
 YEAST CULTURE (CXYST)  
 Zika Virus Nucleic Acid Detection Test, Urine (ZIKVU) –**MUST BE COLLECTED ALONGSIDE PLASMA (EDTA/NAC/ACD-A) TUBE**

**DO NOT USE FOR:**

- Gastrointestinal Pathogens Nucleic Acid Test (GIPC3) – **USE CARY-BLAIR TRANSPORT VIAL**
- OVA & PARASITE EXAM (OP) – **USE PARA PAK FORMALIN AND PVA VIALS**
- STOOL CULTURE (CXSTO, CXARO, CXVIB, CXYER) – **USE CARY-BLAIR TRANSPORT VIAL**

URINALYSIS W/ CULTURE REFLEX, (UACXI) – **MUST COLLECT URINALYSIS AND CULTURE TUBES SEPARATELY**  
 URINE CULTURE (CXURN)

Acid-Fast Bacillus (AFB) Culture and AFB Stain (CXAFB)  
 ANAEROBIC CULTURE (CXANA)  
 Bacterial Culture and Gram Stain (CXBS)  
 Bacterial Culture with Anaerobe and Gram Stain (CXBSN)  
 FUNGAL CULTURE (CXFUN)  
 Mycobacterium tuberculosis/Rifampin Resistance PCR (TBRIF)  
 Sterile Body Fluid Culture with Gram Stain (CXFLD)  
 Sterile Body Fluid Culture with Anaerobe and Gram Stain (CXFLN)



**CAPPED SYRINGE  
 NEEDLE MUST BE  
 REMOVED**



**COPAN eSWAB**

ANAEROBIC CULTURE (CXANA)  
 Bacterial Culture and Gram Stain (CXBS)  
 Bacterial Culture with Anaerobe and Gram Stain (CXBSN)  
 CANDIDA AURIS (CANAT)  
**FUNGAL CULTURE (CXFUN) – STERILE CUP PREFERRED**  
 Group A Streptococcus Culture, Throat (CXGRA)  
 Group B Streptococcus Nucleic Acid Detection Test, Vaginal/Rectal (GBS2)  
 Methicillin Resistant Staphylococcus aureus (MRSA) Screen, Non-nasal (CXMRN)  
 Methicillin Resistant Staphylococcus aureus (MRSA), Surveillance, Nasal (CXMRS)  
 Neisseria gonorrhoeae Culture (CXGON)  
 YEAST CULTURE (CXYST)

**DO NOT USE FOR:**

- CHLAMYDIA/GC PCR (CTNGG, CTGNR, CTNGT) – **USE ROCHE COBAS SWAB KIT**
- FLU, COVID-19, OR VIRAL PCR (HVZPL, FLUCO, FLUPC, RPNA3, COV19) – **USE COPAN UTM SWAB (PINK)**
- M. TUBERCULOSIS RIFAMPIN RESISTANCE PCR (TBRIF) – **USE STERILE CUP, SYRINGE, OR SUCTION TRAP**
- RAPID GROUP A STREP ANTIGEN (RPST) – **USE DRY RAYON OR POLYESTER SWAB**
- STERILE BODY FLUID CULTURE W/GRAM STAIN (CXFLD/CXFLN) – **USE STERILE CUP OR SYRINGE**
- TISSUE CULTURE W/GRAM STAIN (CXTIS/CXTIN) – **USE STERILE CUP**
- VAGINITIS PCR (BCTVP) – **USE HOLOGIC APTIMA KIT**



**HOLOGIC APTIMA KIT**

VAGINITIS PCR (BCTVP)

**DO NOT USE FOR ANY OTHER TEST**



**ROCHE COBAS PCR SWAB**

CHLAMYDIA/GC PCR-GENITAL (CTNGG)  
 CHLAMYDIA/GC PCR-RECTAL (CTNGR)  
 CHLAMYDIA/GC PCR-THROAT (CTNGT)

**DO NOT USE FOR:**

- CHLAMYDIA/GC PCR-URINE – **USE COBAS PCR URINE SAMPLE KIT OR STERILE CUP**
- HPV PCR (HPVC) – **USE THIN PREP VIAL**
- HSV / VZV PCR (HVZPL) – **USE COPAN UTM SWAB**
- VAGINITIS PCR (BCTVP) – **USE HOLOGIC APTIMA KIT**



**SUCTION TRAP**

Acid-Fast Bacillus (AFB) Culture and AFB Stain (CXAFB)  
 ASPERGILLUS GALACTOMANNAN AG BRONCH WASH (ASPGW)  
 Bronchial Washing Culture, Quantitative (CXQBR)  
 FUNGAL CULTURE (CXFUN)  
 Mycobacterium tuberculosis/Rifampin Resistance PCR (TBRIF)  
 Pneumonia Pathogen Nucleic Acid Detection Test (BALPC)  
 RESPIRATORY CULTURE W/GRAM STAIN (CXRES)




**UTM SWAB KIT (RED OR  
 WHITE CAP)**

HSV-1/2 and VZV Nucleic Acid Detection Test, Lesion (HVZPL)  
 Influenza A/B and RSV Nucleic Acid Detection Test (FLUPC)  
 Mpox Virus Nucleic Acid Detection Test, Lesion (MPXNA)  
 PERTUSSIS DNA DETECTION (BORPC)  
 SARS-CoV-2 (COVID-19) Nucleic Acid Detection Test (COV19)  
 UPPER RESPIRATORY PATHOGEN NUCLEIC ACID DETECTION TEST (RPNA3)

**DO NOT USE FOR:**

- ANY BACTERIAL, AFB, OR FUNGAL CULTURE – **USE STERILE CUP, COPAN eSWAB, CAPPED SYRINGE, OR SUCTION TRAP**
- CHLAMYDIA/GC PCR (CTNGG, CTNGR, CTNGT) – **USE ROCHE COBAS PCR SWAB KIT**
- RAPID GROUP A STREP ANTIGEN (RPST) – **USE DRY RAYON OR POLYESTER SWAB**
- VAGINITIS PCR (BCTVP) – **USE HOLOGIC APTIMA KIT**

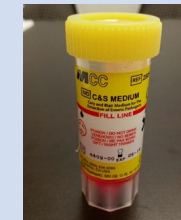

**CARY-BLAIR STOOL  
 TRANSPORT MEDIUM**

AEROMONAS CULTURE (CXARO)  
 VIBRIO CULTURE (CXVIB)  
 YERSINIA CULTURE (CXYER)  
 STOOL CULTURE (CXSTO)  
 Giardia/Cryptosporidium Antigen, Stool (GCEIA)  
**OR**  
 Gastrointestinal Pathogens Nucleic Acid Test (GIPC3)

**DO NOT USE FOR:**

- C. DIFFICILE TOXIN PCR (CDPCI) – **USE STERILE CUP**
- IMMUNOCHEMICAL OCCULT BLOOD (IFOB) – **USE OC-LIGHT S FIT STRIPS**
- STOOL H PYLORI ANTIGEN (HELO) – **USE STERILE CUP**
- WBC SMEAR (WBCSM) – **USE STERILE CUP**



**PARA-PAK  
 PVA & FORMALIN VIALS**

Giardia/Cryptosporidium Antigen, Stool (GCEIA)  
 OVA & PARASITE EXAM, STOOL (OP)

**DO NOT USE FOR:**

- C. DIFFICILE TOXIN PCR (CDPCI) – **USE STERILE CUP**
- Gastrointestinal Pathogens Nucleic Acid Test (GIPC3) – **USE CARY-BLAIR VIAL**
- IMMUNOCHEMICAL OCCULT BLOOD (IFOB) – **USE OC-LIGHT S FIT STRIPS**
- STOOL CULTURE (CXSTO) – **USE CARY-BLAIR VIAL**
- STOOL H PYLORI ANTIGEN (HPYG) – **USE STERILE CUP**
- WBC SMEAR (WBCSM) – **USE STERILE CUP**



**COBAS PCR URINE SAMPLE KIT**

CHLAMYDIA/GC PCR-URINE (CTGCU)

**DO NOT USE FOR ANY OTHER TEST**



**OC-LIGHT FIT TUBES**

IMMUNOCHEMICAL OCCULT BLOOD (IFOB)

**DO NOT USE FOR ANY OTHER TEST**